

TECHNICAL SPECIFICATIONS

Frequency	27,000~32,000 Hz
Rpm	0~32,000
Dimensions	D260 * W140 * H210 mm
Volume	1,000ml
Weight	3.1Kg
E-Type Handpiece Weight	250 g
Piezo Handpiece Weight	45 g
Input	AC 115V/60 Hz or AC 230V/50 Hz



Apixia Corporation Company Address: 7F-1, 128, Sec.2, Dunhua S. Rd., Da-An Dist., Taipei City Taiwan R.O.C.
Office: +886-2-27022608 Fax: +886-2-27554919 Email: vet@apixia.asia www.apixiavet.com


Apixia Inc. Company Address: 18688 San Jose Ave. Industry, CA 91748 USA
Office: 1-626-810-9201 Fax: 1-626-810-9207 Email: info@apixia.com www.apixiavet.com

Apixia Japan Company Address: 1324, Ooida, Daian-cho, Inabe-shi, Mie-ken 511-0274, Japan
Office: 81-594-41-3080 Fax: 81-594-41-3081 Email: sales3@apixia.asia www.apixia.co.jp

BC 160
Rev. D-Aug/06/2021




SCALEX® 890 PLUS ULTRASONIC PIEZO SCALER and MICROMOTOR COMBINATION

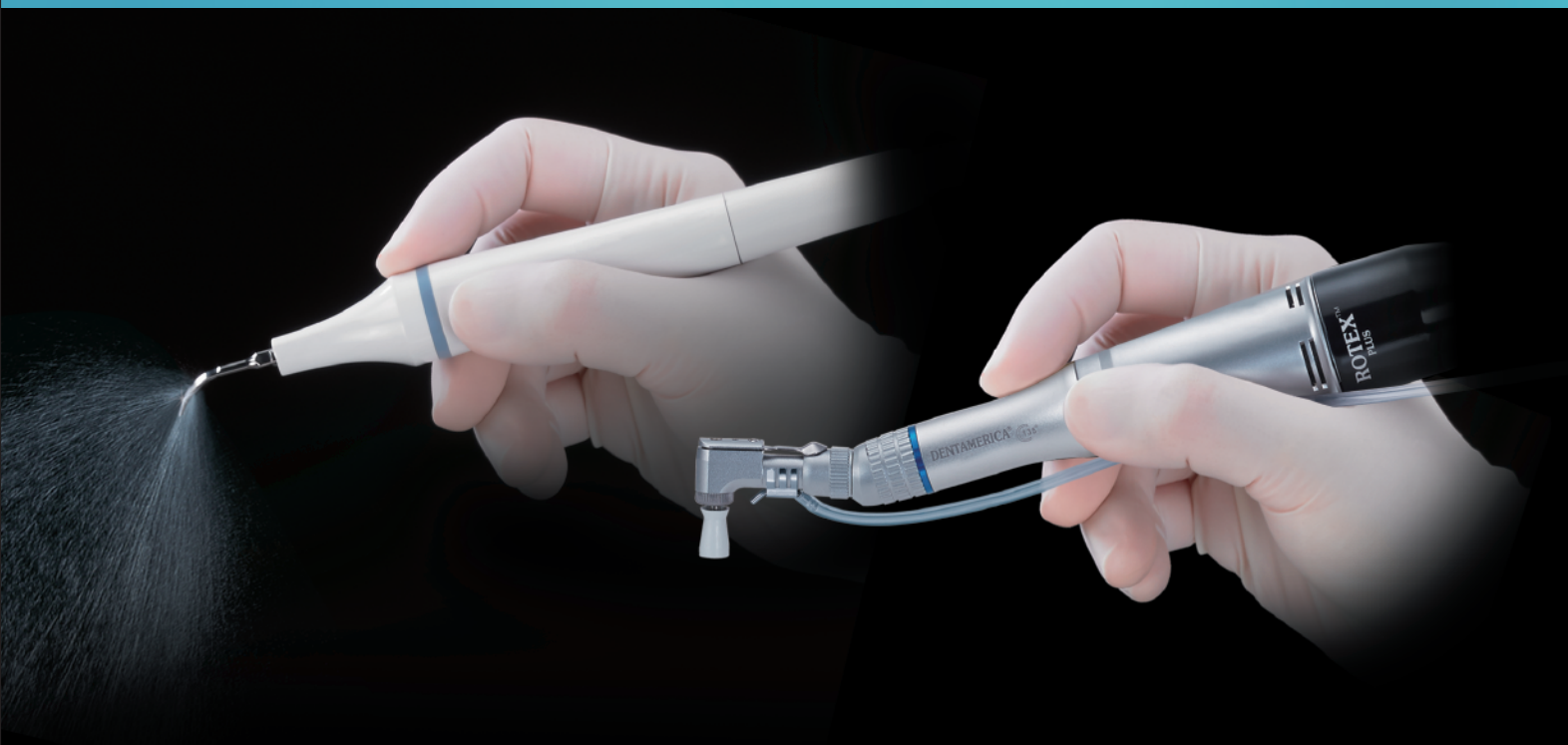

SOLUTION TANK

- Large 1,000 ml solution tank with internal pump allows for convenient portability.
- Any type of irrigant can be used during procedures to meet all clinical needs.

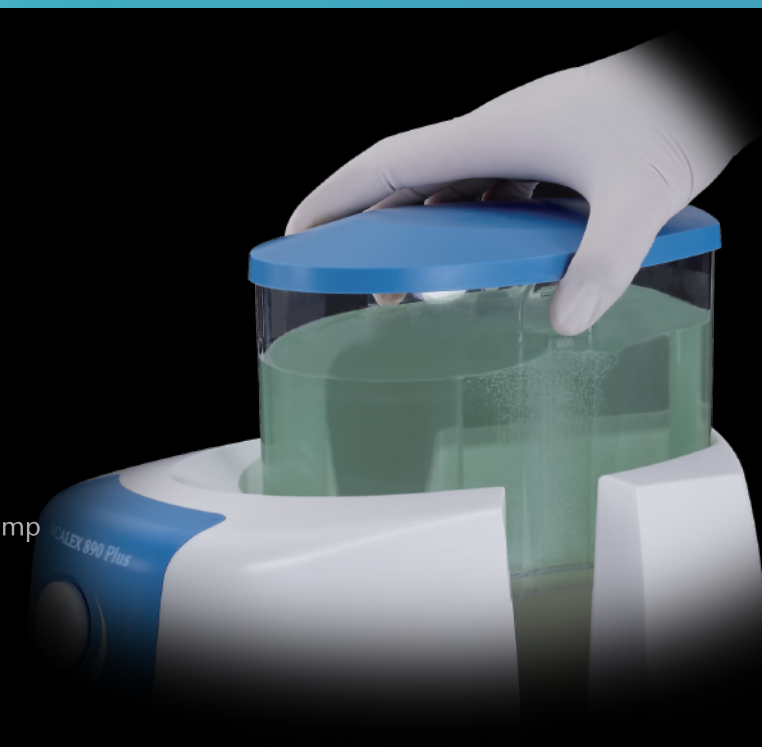

TWO ANIMALS IN ONE FOR ALL OF YOUR PET HYGIENE NEEDS





- ANIMAL HELPER:** Prevents periodontitis, tooth loss and bad breath in pets.
- MULTI-FUNCTION:** Unique and versatile dual purpose veterinary device.
- DUAL EFFECTIVENESS:** Contains a self-contained piezo ultrasonic scaler and an electric micro motor in one compact unit.
- SINGLE TOUCH:** Conveniently switch between scaler and micromotor at the touch of a button.
- EASY ACCESS:** Universal foot control accessible from 360°.

Self-contained Scaler

- Lightweight and streamlined autoclavable piezo handpiece.
- A variety of piezo scaler tips are available.

Electric Micromotor

- Straight and contra angle attachments included for prophylaxis and a variety of other functions.
- Variable speed ranges from 0 to 32,000 RPM.
- Automatic cooling fan for continuous high speed performance.

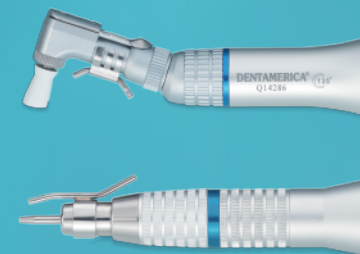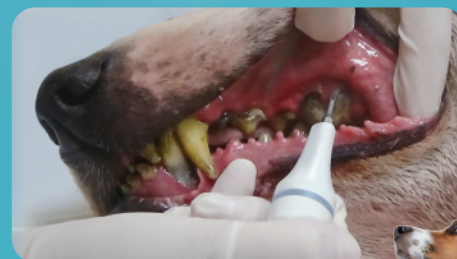

STRAIGHT AND CONTRA ANGLE ATTACHMENTS

Application for prophylaxis, grinding, polishing of animal's nails and shell with sterilizable coolant delivery tube.

A VARIETY OF PIEZO SCALER TIPS ARE AVAILABLE.

- Endo Chuck: Root canal preparation and cleaning.
- Universal A: Removal of all supragingival calculus.
- Universal B : Removal supragingival calculus and suitable for interdental cleaning.
- Chisel: Removal supragingival heavy calculus.
- Perio: Removal of subgingival deposits and deep perio calculus.

